Anleitung

Pedalo® PhysioFlip® Art.Nr. 285 9402

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Tel. +49 (0) 7381 - 9357 0
Fax +49 (0) 7381 - 9357 40

info@pedalo.de
www.pedalo.de

Before using the device for the first time
(or before starting with assembly),
please **read instruction completely** and keep it.

consists of:

- Base with oval base plate
- Rubber webbing (628 181)
- Adjustable handle
- Level compensation for support leg

Important information and warning hints:

Warning! The packing is not part of the product and have to be disposed.
Keep plastic bags and films out of the reach of babies and young children – RISK OF SUFFOCATION!

Warning! Use the device exclusively for the purposes described in this user manual.
Even if it is, among other things, a play device for children, the assembly must always be carried out or supervised by an adult.

Warning! Only use device on level ground with even and firm surfaces.
After longer or frequent use, check that the screws are tight and retighten if necessary.
Defective and worn parts should be replaced by original replacement parts immediately.

The product is made of natural beech wood. Differences in color, structure, grain, small included knots or scars, prove the genuineness of the material and therefore show no reduction in quality.


Maintenance and cleaning:
Can be wiped off with a damp cloth. Do not use solvent-based, acidic or abrasive cleaners.
Protect item from wet.

If required, the cork layer can be rubbed with a disinfecting cloth (please observe the respective manufacturer's instructions).


This device is meant for using from persons at the age of 12 up.

Maximum weight capacity 120 kg

Assembly elastic webbing

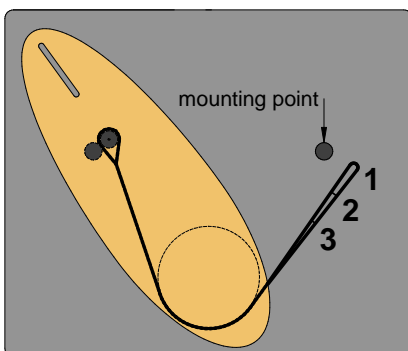
The Pedalo-Physio Flip® can be used left or right. Insert the rubber webbing as shown in the illustration and fasten it to the respective attachment point. Please note that the webbing does not protrude at the attachment points, this could lead to annoying noises and wear of the webbing during operation. When changing direction, the webbing is pushed to the other side under the footrest and fixed again. Three different preloads can be selected
1 = weak preload
2 = medium preload
3 = high preload

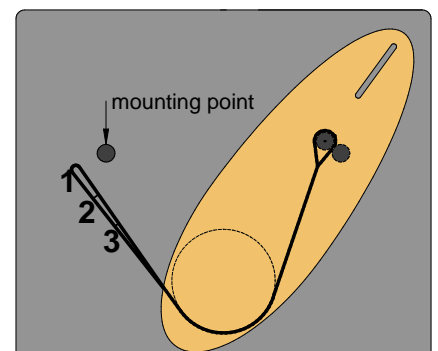Assembly handle

The handle is fixed in the long slot of the oval base plate.
Please refer to the drawing on the right for the assembly sequence of the handle.
The long slot allows the handle to be infinitely adjusted to the distance to the centre of rotation. After adjustment, tighten the handle firmly with the wing nut.

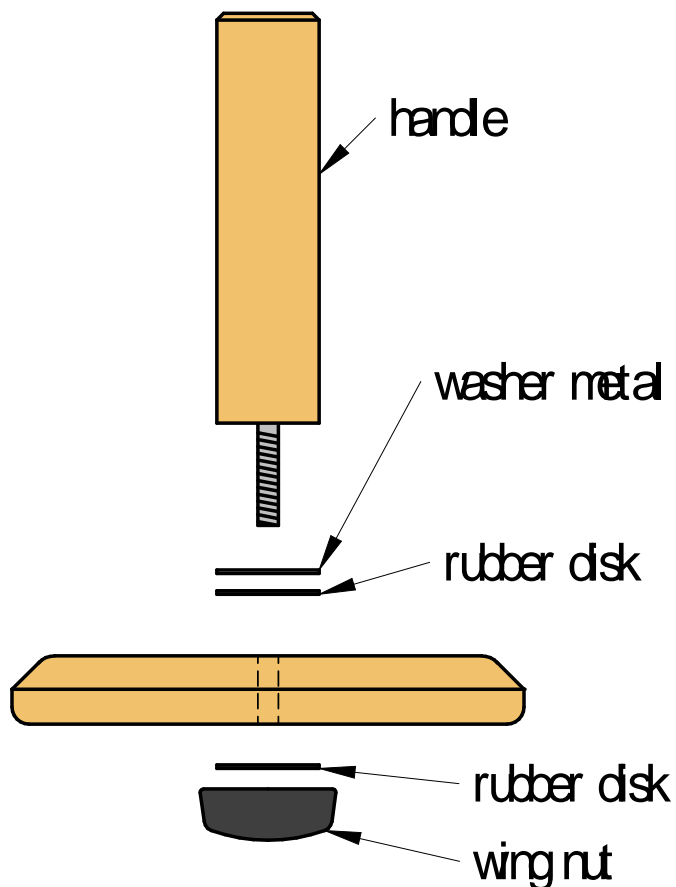Level compensation:

The level compensation is placed close to the base plate for foot exercises while standing or sitting. The support leg can be positioned on it.

Spare parts:

Rubber webbing for the Pedalo-PhysioFlip® (item no. 628 181)

To replace, unhook the rubber strap, then unscrew one of the stop buffers. Simply turn the foot rest until it is possible to reach the second attachment point of the rubber strap from below. Now change the rubber webbing, then swivel the foot rest back into its normal position and screw on the stop buffer again. Finally tighten the rubber strap again and attach it to one of the stops.