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NOTE:


Despite this instruction the enclosed manual needs to be observed carefully. Especially the chapter about safety and dangers, the intended use and the foreseeable misuse needs to be observed.

1 OVERVIEW


- (1) Foot rest (application part type BF)
- (2) Crank disc
- (3) Control and display unit with 2.7" colour screen
- (4) Holder for control and display unit with 2.7" colour screen
- (5) Helix cable
- (6) Therapy grip (application part type BF)
- (7) Height-adjustable foot
- (8) Base unit

2 INSTALLING FOOT RESTS

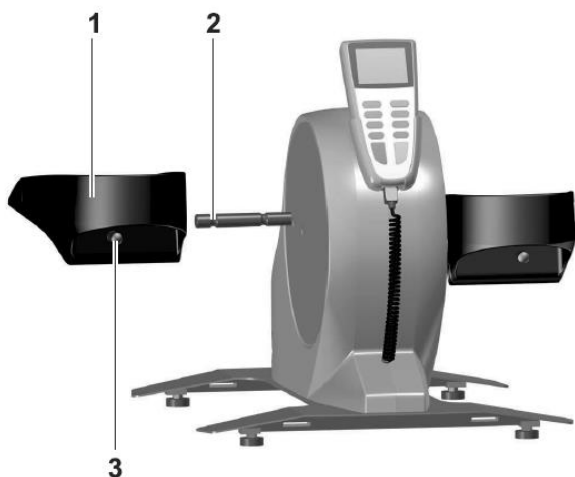
- ▶ Press locking bolt **3**.
- ▶ Slide foot rests **1** onto crank axle **2** up to the stop.
- ▶ Release locking bolt **3**.
- ▶ Make sure that locking bolts audibly click into place.

3 INSTALLING THERAPY GRIPS

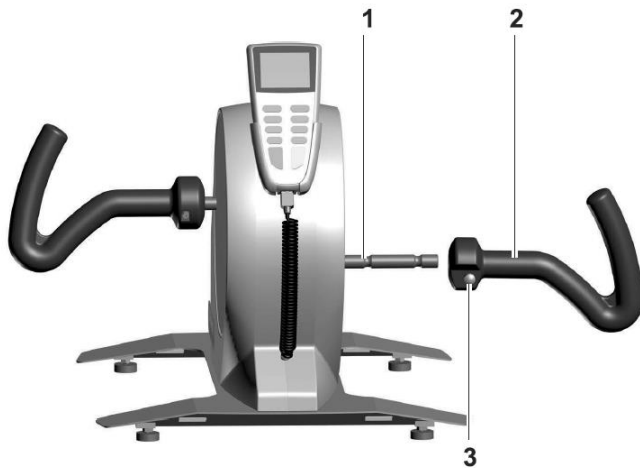- ▶ Press locking bolt **3**.
- ▶ Slide therapy grips **2** onto crank axle **1** up to the stop.
- ▶ Release locking bolt **3**.
- ▶ Make sure that locking bolts audibly click into place.


4 ADJUSTING HEIGHT-ADJUSTABLE FEET

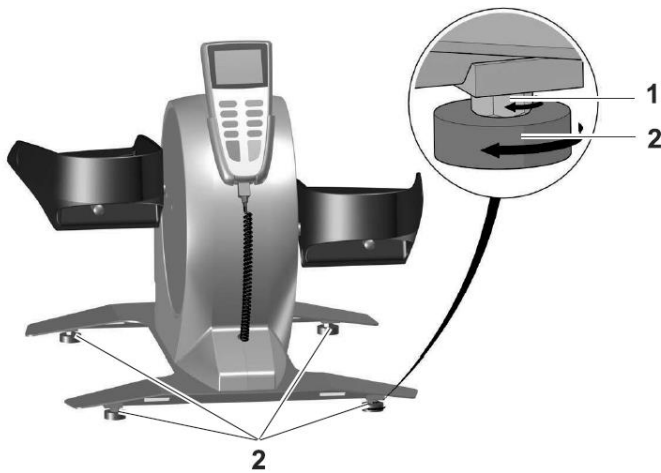To level out unevenness of the floor, it is possible to adjust the height of the device's feet to improve the stability and slip resistance of the THERA-Trainer mobi.


Adjust feet (4) as follows:

- ▶ Make sure that the THERA-Trainer mobi is without power, pull power plug if necessary.
- ▶ Unscrew lock nuts **1** using a suitable tool.
- ▶ Adjust feet **2** to required height.
- ▶ Make sure that the THERA-Trainer mobi is levelled.
- ▶ Tighten lock nuts and verify tightness.
- ▶ Make sure that feet of device are dry and free of dust.
- ▶ The manufacturer recommends cleaning the device's feet before every training session.

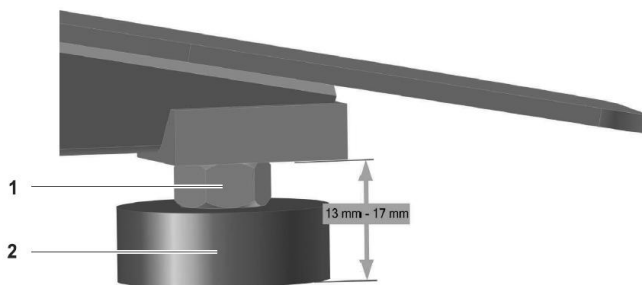The height-adjustable feet serve exclusively to level out unevenness. The maximum adjusting height between surface and metal frame is 13 mm to 17 mm. If the adjusting range of 4 mm max. is not respected, it is possible for the lock nut to damage the thread.

5 CONNECTING THE CONTROL AND DISPLAY UNIT WITH 2.7" COLOUR SCREEN

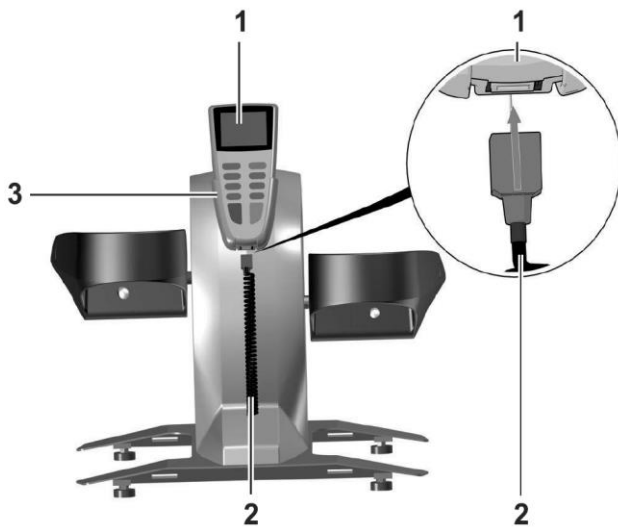
- ▶ Make sure that control and display unit is without power.
- ▶ Plug helix cable 2 into control and display unit 1.

Control and display unit is correctly connected when plug connector audibly clicks into place.

- ▶ Place control and display unit in the intended holder 3.


6 CONNECTING THE POWER CABLE

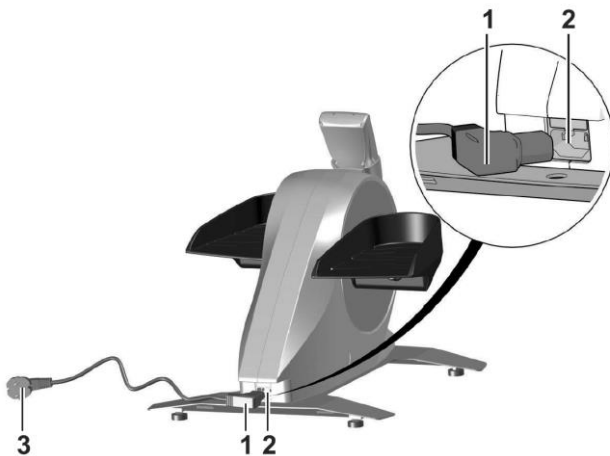- ▶ Use exclusively the power cable supplied by the manufacturer.
- ▶ Connect plug 1 to port 2.
- ▶ Connect power plug 3 to socket outlet.

The THERA-Trainer mobi is ready for use.

7 CONTROL AND DISPLAY UNIT WITH 2.7" COLOUR SCREEN

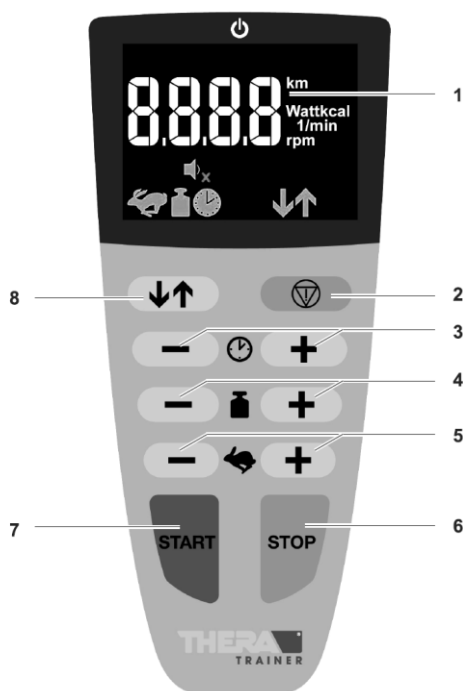- (1) Training parameters display
- (2) EMERGENCY STOP button
- (3) Training time buttons
- (4) Training resistance buttons
- (5) Number of revolutions buttons
- (6) STOP button
- (7) START button
- (8) Direction of rotation button

Safe handling of the control and display unit:

- ▶ Make sure that you can reach the control and display unit at any time to stop THERA-Trainer mobi in emergency situations.
- ▶ Make sure that helix cable of the control and display unit with 2.7" colour screen does not get tangled in the therapy grips or foot rests.

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medica Medizintechnik GmbH