
PROFESSIONAL CARE









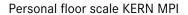

Personal scale with particularly large, flat weighing plate to make it easier to step on and give maximum safety

Features

- This personal scale is the ideal measuring instrument for elderly care, paramedics or private sector, and can weigh overweight people up to 200 kg
- Very large weighing plate
- Easy and hygienic cleaning
- In the scale can be easily transported using a handle and does not require much storage space
- 2 Particularly flat design facilitates access
- Server Very large platform with non-slip and wear-resistant plastic-surface
- Easy to operate: simply tap on it stand on it – read the weight
- Integrated AUTO-OFF function to preserve the batteries
- · Weight displayed in kg or lb

STANDARD

UNIT

BATT 1 DAY

- · Large display that is easy to read
- Delivered in single design packaging
- Note: The models are only delivered in a set of 5 units. That means the price refers to a delivery of 5 units. Cannot be delivered individually

Technical data

- Large LCD display, digit height 30 mm
- Overall dimensions W×D×H
- Batteries included, 4×1,5 V AAA, operating time up to 20 \mbox{h}
- Net weight approx. 8 kg

Accessories

 Cleaning cloths, alcohol-free cloths for disinfectant wiping, quick acting, based on modern quaternary ammonium compounds and effective against papova viruses. Particularly gentle on materials, and very well suited for disinfecting products which are sensitive to alcohol. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 100 pcs., size 20×22 cm per cloth, KERN MYC-01

Model	Weighing capacity	Readability	
KERN	[Max] kg	[d] kg	
MPI 200K-1S05	200	0,1	

KERN & SOHN GmbH · Ziegelei 1 · 72336 Balingen · Germany · Tel. +49 7433 9933 - 0 · www.kern-sohn.com · info@kern-sohn.com

MEDICAL SCALES 2023

KERN PICTOGRAMS



• 6888 •

```
Adjusting program CAL:
          For quick setting up of the balance's accuracy. External adjusting weight
          required
           Memory:
           Balance memory capacity, e.g.
MEMORY
          for article data, weighing data, tare weights, PLU etc.
           Data interface RS-232:
RS 232
```

To connect the balance to a printer, PC or network



RS-485 data interface:

To connect the balance to a printer, PC or other peripherals. Suitable for data transfer over large distances. Network in bus topology is possible



Bluetooth* data interface:

To transfer data from the balance to a printer, PC or other peripherals



WIFI data interface: To transfer data from the balance to a printer, PC or other peripherals



Control outputs (optocoupler, digital I/O): To connect relays, signal lamps, valves, etc.



Statistics:

using the saved values, the device calculates statistical data, such as average value, standard deviation etc.



PC Software: to transfer the measurements from the device to a PC



GLP/ISO log: With date and time. Only with KERN printers



KERN Communication Protocol (KCP):

It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems



Piece counting: Reference quantities selectable. Display can be switched from piece to weight



З

Totalising level A: The weights of similar items can be added together and the total can be printed out

Weighing units: Can be switched to e.g. nonmetric units. Please refer to website for more details




Weighing with tolerance range


(Check weighing) Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model




ZERO: Resets the display to "0"

Hold function:

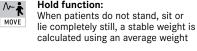

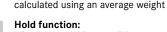

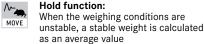



Protection against dust and water splashes IPxx:

The type of protection is shown in the pictogram cf. DIN EN 60529:2000-09, IEC 60529:1989+A1:1999+A2:2013



Suspended weighing: Load support with hook on the underside of the balance



Battery operation: Ready for battery operation. The battery type is specified for each device



Rechargeable battery pack: Rechargeable set



Battery operation rechargeable Prepared for a rechargeable battery operation

Universal plug-in power supply: with universal input and optional

230V/50Hz in standard version for EU.

On request GB, AUS or USA version

input socket adapters for

Plug-in power supply:

A) EU, CH B) EU, CH, GB, USA

available



230 V



Integrated power supply unit: Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, AUS or USA on request



Weighing principle: Strain gauges Electrical resistor on an elastic deforming body



Peak hold function: capturing a peak value within a measuring process



Push and Pull: the measuring device can capture tension and compression forces



KERN & SOHN GmbH · Ziegelei 1 · 72336 Balingen · Germany · Tel. +49 7433 9933 - 0 · www.kern-sohn.com · info@kern-sohn.com



360° rotatable microscope head



Monocular Microscope: For the inspection with one eye

KERM



Binocular Microscope: For the inspection with both eyes

Trinocular Microscope:

For the inspection with both eyes and the additional option for the connection of a camera



Abbe Condenser: With high numerical aperture for the concentration and the focusing of light



Halogen illumination: For pictures bright and rich in contrast



LED illumination:

Cold, energy-saving and especially long-life illumination



Fluorescence illumination for compound microscopes: With 100 W mercury lamp and filter



Fluorescence illumination for compound microscopes: With 3W LED illumination and filter



PH

Phase contrast unit: For a higher contrast



Darkfield condenser/unit: For a higher contrast due to indirect illumination

∰ POLAR

Polarising unit: To polarise the light



Infinity system: Infinity corrected optical system



Automatic temperature compesation: For measurements between 10 °C and 30 °C



Verification possible: The time required for verification is specified in the pictogram

1 DAY

Package shipment:

Pallet shipment:

pictogram

The time required for internal shipping preparations is shown in days in the pictogram

The time required for internal shipping

preparations is shown in days in the



*The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by KERN & SOHN GmbH is under license. Other trademarks and trade names are those of their respective owners.